



MICRON[®]

MICRON Animal Stand - Pro

On the web: <https://phoenixmicron.com>

Email us at: sales@phoenixmicron.com



MICRON Animal Stand - Pro - Mouse

With the **MICRON Animal Stand - Pro**, for the first time, inhalation anesthetic delivery is integrated into a platform specifically designed for *in vivo* ophthalmic imaging of small animals.



MICRON Animal Stand - Pro - Rat

Supply and venting connections are located on the base of the animal stand, away from the animal's eyes and the objective lens for simple connection to your delivery system

Fully adjustable platform with:

- Precise height micro-adjustments
- Horizontal plane rotation
- Adjustable horizontal tilt functionality






MICRON Animal Stand - Pro

On the web: <https://phoenixmicron.com> Email us at: sales@phoenixmicron.com



An innovative mask design allows for optimal presentation of the subject's left and right eyes to the camera.

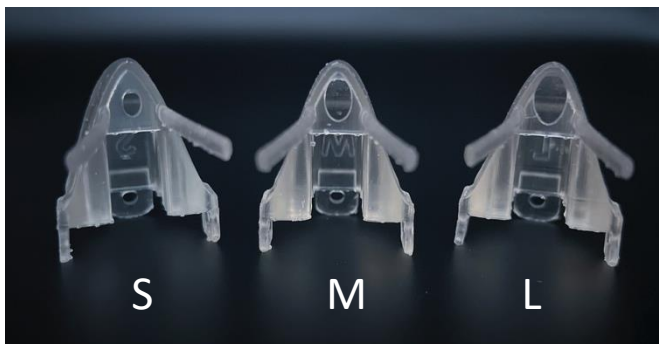

The **WhiskerWipes® Animal Mask** is available in separate sizes for mice and rats and snaps securely on the animal platform. This provides integrated inhalation anesthesia delivery and venting, enabling sustained animal sedation.

WhiskerWipes® (patent pending) arms hold whiskers away from the objective lens, simplifying the imaging process and minimizing artifacts.

The soft, snug-fitting masks are made of non-allergenic, biocompatible materials that are both washable and can be easily replaced to comply with laboratory infection control requirements.





MICRON[®]

MICRON Animal Stand - Pro


On the web: <https://phoenixmicron.com> Email us at: sales@phoenixmicron.com

MICRON Animal Stand – Pro with MICRON Heater



Maintain optimum subject animal body temperature for long sessions by adding the fully integrated **Heater for the MICRON Animal Stand – Pro**.

- The easily-cleaned aluminum animal platform becomes the heater, eliminating the need for an additional component placed on top of the animal stand.
- Wiring is routed through the inside of the animal stand to keep the space around the animal eyes and the MICRON camera free of wire clutter.

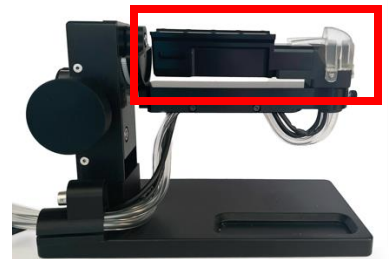

Tubing connection size:

- 1/4" Outer Diameter
- 1/8" Inner Diameter





MICRON Animal Stand - Pro

On the web: <https://phoenixmicron.com> Email us at: sales@phoenixmicron.com

Upgrade to the MICRON Animal Stand - Pro



MICRON Mouse - Pro Upgrade Kit



MICRON Rat - Pro Upgrade Kit

Upgrade Kits are available for the MICRON Animal Stand-Standard!

Upgrade your existing MICRON Animal Stand – Standard to a MICRON Animal Stand – Pro for MOUSE or for RAT.

- User-installable Pro upgrade for MICRON Animal Stand - Standard
- Integrated inhalant anesthesia delivery system
- Includes 3 WhiskerWipes® Animal Masks and installation tools
- Compatible with the Heater for MICRON Animal Stand – Pro
- Platforms options for MICE or RATS
- Tubing connection size:
 - 1/4" Outer Diameter
 - 1/8" Inner Diameter





MICRON Animal Stand - Pro

On the web: <https://phoenixmicron.com> Email us at: sales@phoenixmicron.com



Compatible with MICRON 5 and MICRON IV Cameras
as well as older MICRON III systems

1000+

Published papers that
incorporate MICRON
data

15+

Years of experience
innovating patented
small animal imaging
technology

1

Integrated multi-
modality **system**

8

Imaging modalities,
designed for exacting
small animal ophthalmic
research

